


Bante 920E Benchtop pH Meter

USER MANUAL

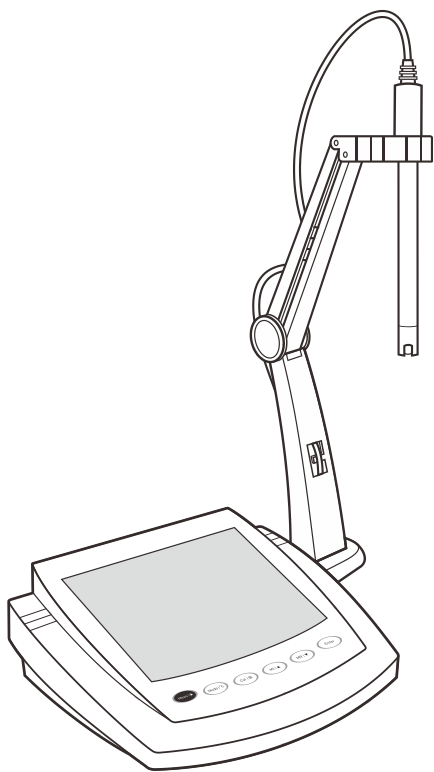

Introduction

Thank you for selecting the 920E benchtop pH meter. This user manual provides a step-by-step guide to help you operate the meter, please carefully read the following instructions before use. Any use outside of these instructions may invalidate your warranty and cause permanent damage to the meter.

Environmental Conditions

Before unpacking, ensure that current environmental conditions meet the following requirements.

- Relative humidity is less than 80%
- Ambient temperature between 0°C (32°F) and 50°C (122°F)
- No potential electromagnetic interference
- No corrosive gas exists

Packing List

The following list describes all components of the meter. If any items are missing or damaged, contact the supplier immediately.

920E meter

Electrode arm

USB cable


DC 5V power adapter

pH electrode

TP-10K temperature probe

pH buffer reagents 4.01, 7.00, 10.01

Meter Overview

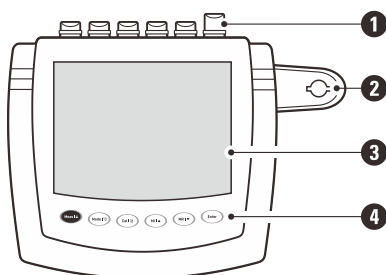


1 Sensor connections

2 Base plate of electrode arm

3 Display

4 Membrane keypad

Connectors

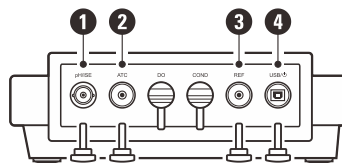


1 Socket for pH or ORP electrode (BNC)

2 Socket for temperature probe (3.5 mm jack)

3 Socket for reference electrode (3.5 mm jack)



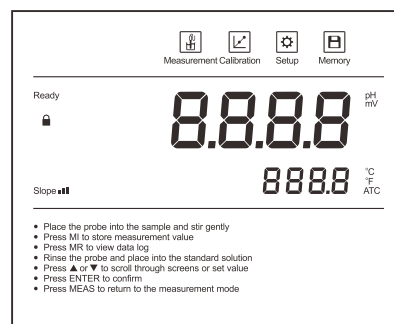



4 USB-B interface to the power adapter or computer

Display




Icon	Description
	Indicates that the meter is in the measurement mode
	Indicates that the meter is in the calibration mode
	Indicates that the meter is in the setup mode
	Indicates that you are viewing the stored readings or a reading is stored into the memory
Ready	Shown when the measurement is stable
	Shown when the reading is locked
Slope ■■■	If the electrode slope exceeds the allowed range after calibration, the icon automatically disappears
ATC	Indicates that the automatic temperature compensation is enabled

Keypad

Key	Function
	<ul style="list-style-type: none"> Switch the meter on or off Lock or unlock the measurement Exit the calibration, settings, data logs and return to the measurement mode
	<ul style="list-style-type: none"> Toggle between the pH and mV modes Press and hold the key to enter the temperature setting
	<ul style="list-style-type: none"> Start calibration Press and hold the key to enter the setup menu
	<ul style="list-style-type: none"> Store current reading to memory Increase value or scroll up through a list of options
	<ul style="list-style-type: none"> View the data log or calibration log Decrease value or scroll down through a list of options
	<ul style="list-style-type: none"> Confirm the calibration or displayed option Press and hold the key to switch the backlight on or off

Installing the Electrode Holder


Take out the electrode arm from the accessory box. The base plate of electrode arm has a circular hole, the electrode arm has a connecting rod. Insert the connecting rod into the circular hole and swivel the electrode arm 90 degrees. The electrode holder is now ready to swing into desired position.

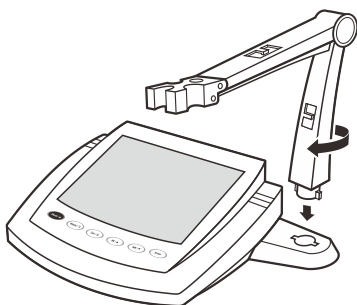


Adjusting the Electrode Arm

After installation, if the electrode arm automatically rises or falls, you are able to adjust the screw until arm locate at any position.

1. Remove the plastic cover from the right side of the electrode arm.

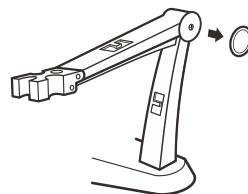

2. Use the screwdriver to tighten the screw moderately.

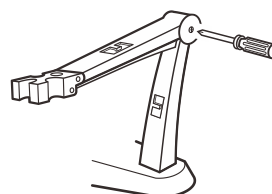


3. Insert the plastic cover to previous position.


Connection

Connecting the pH Electrode

- 1.1 Take out the electrode from packaging. Follow the steps below to place electrode into the left or right side of the electrode arm.

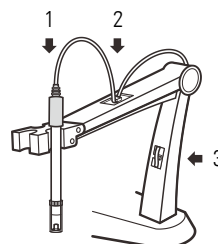


- 1.2 Insert BNC connector into the connector socket labeled pH/ISE. Rotate and push the connector clockwise until it locks.

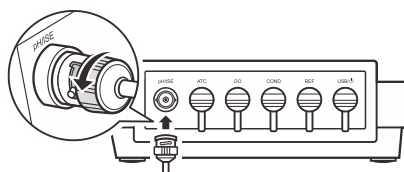


- 1.3 After the connection is completed, DO NOT pull on the cable. Always make sure that the connector is clean and dry.

Connecting the Temperature Probe

- 2.1 Place the temperature probe into the circular hole located at the center of the electrode arm.

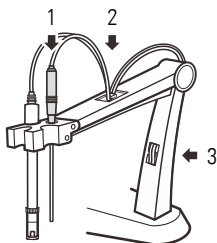


- 2.2 Insert the jack plug to the connector socket labeled ATC. Ensure the connector is fully seated.

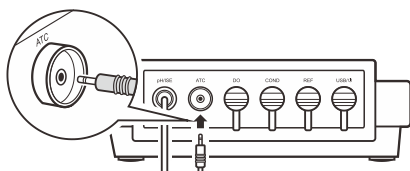


Connecting the Power Adapter

- 3.1 Connect the USB cable to the meter and power adapter.
3.2 Plug the DC 5V power adapter into the wall outlet.

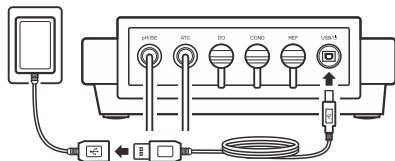


Prior to Use

Remove the protective cap from the bottom of the pH electrode. If tiny air bubbles are present inside the pH-sensitive glass membrane, gently shake the electrode downward to remove air bubbles.

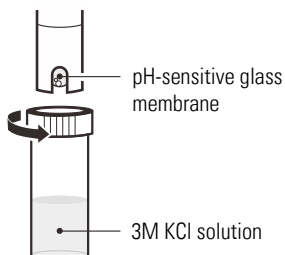


Preparation of pH Buffer Solutions

1. The meter is packaged with pH 4.01, 7.00, 10.01 buffer reagents required for calibration. Half fill a 250 ml volumetric flask with distilled water and add the pH 7.00 buffer reagent.

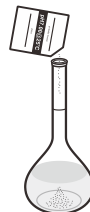

2. Swirl the volumetric flask gently to dissolve the reagent and fill to the mark with distilled water.

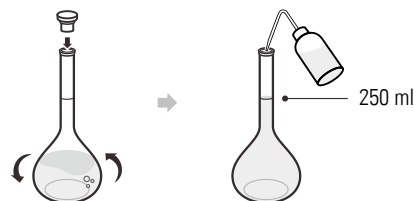


3. Cap and upend the volumetric flask several times to mix the solution.



- Preparation of pH 4.01 and 10.01 buffer solutions are the same as above.
- Prepared buffer solution should be stored in hermetically sealed glass container and avoid direct sunlight.

Switching the Meter On and Off

- Press the **Meas** key for about 5 seconds and release to switch on the meter.
- Press and hold the **Meas** key to switch off the meter.

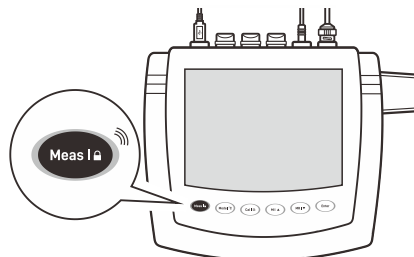


Meter Setup

The 920E meter contains an integrated setup menu for customizing the displayed option to meet measurement requirements. The following table describes the functions of each menu item.


Menu Item	Option and Description	
CAL	Calibration Points Set the number of calibration points.	
	1	1 point
	2	2 points
	3	3 points (default)
BUF	pH Buffer Group Set the pH buffer group for calibration and auto-recognition.	
	USA	USA (default)
	NIST	NIST
UNIT	Measurement Unit Set the default temperature unit.	
	°C	Degrees Celsius (default)
	°F	Degrees Fahrenheit
HOLD	Auto-Hold If enabled, the meter will automatically sense and lock the measurement endpoint.	
	YES	Enable
	NO	Disable (default)
OFF	Auto-Power Off If enabled, the meter will automatically switch off if no key is pressed within 30 minutes.	
	YES	Enable
	NO	Disable (default)
CLR	Clear Stored Data Delete all data logs in the memory.	
	YES	Enable
	NO	Disable (default)
RS	Factory Reset Reset the meter to factory default settings. Note, the meter must be recalibrated.	
	YES	Enable
	NO	Disable (default)

Setting a Default Option


1. In the measurement mode, press and hold the **Cal I** key to enter the setup menu.
2. Press the **▲** / **▼** key to select a menu item.

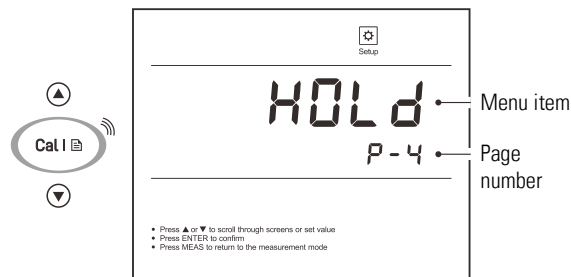

3. Press the **Enter** key, the meter shows the current option.

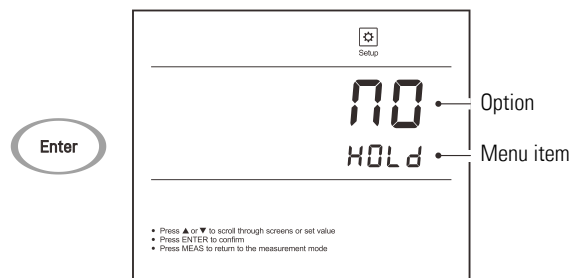

4. Press the **▲** / **▼** key to select a desired option.
5. Press the **Enter** key to save and return to the measurement mode.

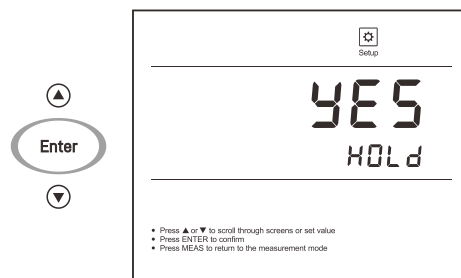

To exit the setup menu without saving changes, press the **Meas** key.

Temperature Compensation

For better accuracy, we recommend the use of either a sensor with a built-in or a separate temperature probe. The meter will calculate the pH slope with measured temperature and show the temperature compensated readings.

Automatic Temperature Compensation

Connect the temperature probe to meter, the ATC icon appears on the display, the meter is now switched to the automatic temperature compensation mode.



Refer to the *Connecting the Temperature Probe* section on page 3.

Manual Temperature Compensation

If the meter does not detect a temperature probe, the degrees Celsius icon (°C) will show on the display indicating the meter is switched to the manual temperature compensation mode. To set the temperature value follow the steps below.

1. Press and hold the **°C** key to enter the temperature setting.
2. Press the **▲** / **▼** key to modify the temperature value.
3. Press the **Enter** key to save.



Press and hold the **▲** / **▼** key will make the value change faster.

pH Calibration

The 920E meter allows up to 3 points pH calibration. We recommend that you perform at least 2 points calibration for high accuracy measurement. The meter will automatically recognize and calibrate to following standard buffer values.


USA Standard Buffers	pH 4.01, 7.00, 10.01
NIST Standard Buffers	pH 4.01, 6.86, 9.18

Single point calibration should only be carried out with the pH 7.00 or 6.86, otherwise calibration will not be accepted.


Make sure to calibrate the meter when attaching a new pH electrode or during first use. Do not reuse the buffer solutions after calibration, contaminants in solution will affect the calibration and eventually the accuracy of the measurement.

Single Point Calibration


- 1.1 Ensure that the meter is in the pH measurement mode and you have selected 1 point calibration in the setup menu.
- 1.2 Press the **Cal** key, the meter shows 7.00/CAL 1 or 6.86/CAL 1, depending on the selected pH buffer group.

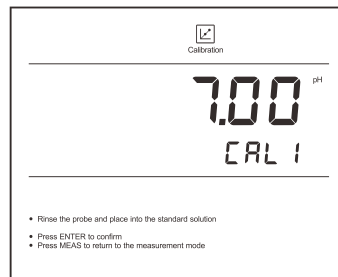

- 1.3 Rinse the pH electrode with distilled water, place the electrode (and temperature probe) into the pH 7.00 buffer solution, stir gently to create a homogeneous solution.

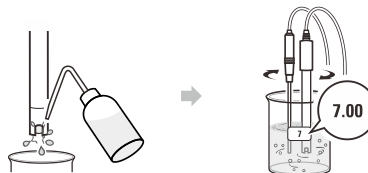


- 1.4 Press the **Enter** key, the Calibration icon begins flashing.

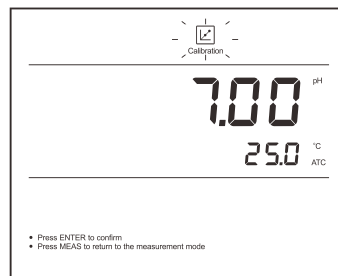


- 1.5 When the reading has stabilized, the meter will show *End* and return to the measurement mode.

Multipoint Calibration


- 2.1 Ensure that you have selected 2 or 3 points calibration in the setup menu.
- 2.2 Repeat steps 1.2 through 1.4 above. When the first calibration point is completed, the display will show ----/CAL 2. The meter prompts you to continue with second point calibration.

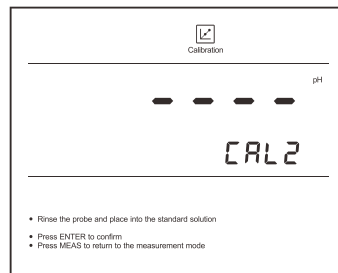


- 2.3 Rinse the pH electrode with distilled water, place the electrode (and temperature probe) into the next buffer solution (e.g., pH 4.01).

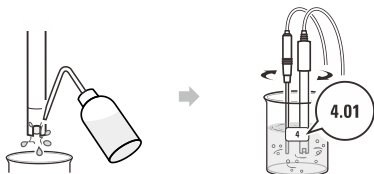

The meter will automatically recognize the buffer solution and begin the calibration, the Calibration icon continuously flashing.

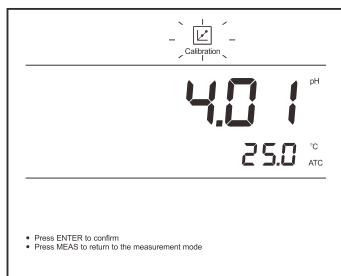


- 2.4 When the reading has stabilized, the display will show ----/CAL3. The meter prompts you to continue with third point calibration.
- 2.5 Repeat the step 2.3 above until the meter shows **End**. Calibration is completed.

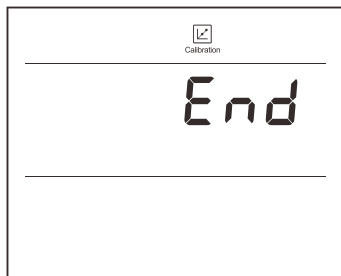


- During the calibration, if the display shows ---- indicating the meter is waiting for recognizing the pH buffer solution.
- If the display shows **ERR** indicating the measured mV value for the current calibration point deviates by more than 60 mV from the theoretical value of the pH buffer. The calibration will not be accepted. Please check the pH electrode and ensure the buffer solutions are fresh and uncontaminated.
- If the calculated electrode slope is not between 70% to 110%, Slope **■■** icon will disappear from the display. The pH electrode may need to be replaced.
- To exit the calibration without saving changes, press the **Meas** key.

Viewing the Calibration Log


- 3.1 Press the **MR** key in the measurement mode and press the ▼ key until the meter shows **EL E / P - 2** (Electrode/Page 2).
- 3.2 Press the **Enter** key, the meter shows the zero-point offset.

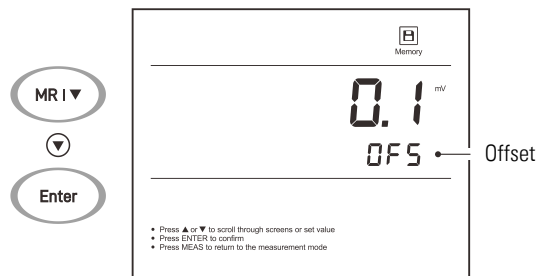

- 3.3 Press the ▼ key to view the calibration point 1
- 3.4 Press the ▼ key to view the calibration point 2.

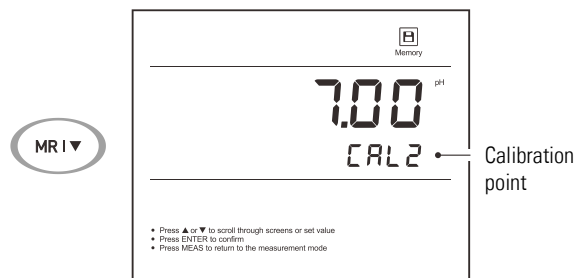

- 3.5 Press the ▼ key to view the electrode slope.

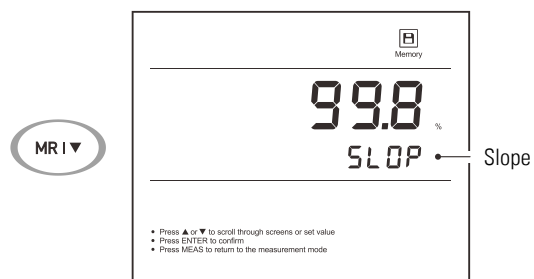

- 3.6 To exit the calibration log, press the **Meas** key.




If the meter is not calibrated, the display will show ---- only.

Temperature Calibration

The 920E meter is supplied with a temperature probe for measurement and temperature compensation. If the measured temperature reading differs from that of an accurate thermometer, the probe needs to be calibrated.


1. Connect the temperature probe to the meter and place into a solution with a known accurate temperature.
2. Press and hold the **°C** key to enter the temperature setting.
3. Press the **▲** / **▼** key to modify the temperature value.
4. Press the **Enter** key to save.

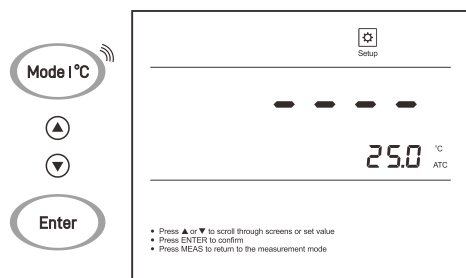


Measurements

pH Measurement

- 1.1 Rinse the pH electrode with distilled water. Place the electrode (and temperature probe) into the sample solution and stir gently. Note, the pH-sensitive glass membrane and liquid junction must be completely immersed into the solution.

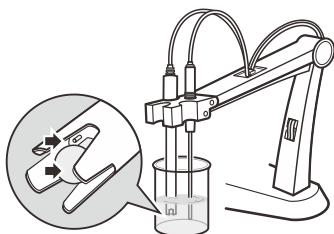

- 1.2 If the Auto-Hold option in the setup menu is enabled, the meter will automatically sense a stable reading and lock measurement, the Ready/🔒 icon appears on the display. Press the **Meas** key to resume measuring. If the option is disabled, the meter will continuously measure and update the readings.
- 1.3 Wait for the measurement to stabilize and record the reading.
- 1.4 When all of the samples have been measured, rinse the electrode according to the instructions in the **Electrode Maintenance**.

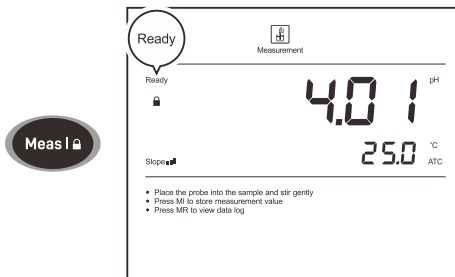

- During the measurement process, never wipe the pH-sensitive membrane as this will cause static interference, blot dry with a lint-free tissue to remove waterdrops on electrode.
- If the meter shows ---- indicating the measurement exceeds the range, remove the electrode from the sample immediately.
- If your sample is pure water, low ionic or low conductivity water, we recommend measuring the pH in the smallest sample volume possible or adding 0.3 ml of the 3M KCl to 100 ml of the sample solution. Note, only high purity KCl can be used.


mV Measurement

- 2.1 Press the **Mode** key unit the display shows measurement unit mV.
- 2.2 Rinse the electrode with distilled water. Place the electrode into the sample solution and stir gently. Record the reading when the measurement is stable.

Data Management


Storing a Measurement Result

In the measurement process, press the **MI** key to store the reading into the memory, the Memory icon appears on the display.

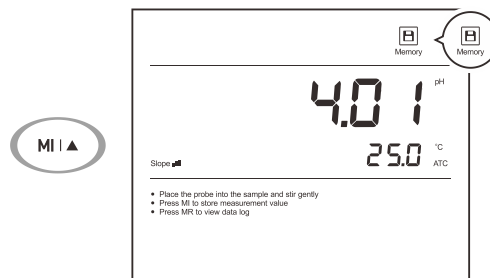


Viewing the Data Logs

- 1.1 Press the **MR** key in the measurement mode, the meter shows **LOG / P - 1** (Log/Page 1).
- 1.2 Press the **Enter** key, the meter shows the serial number of the stored data.

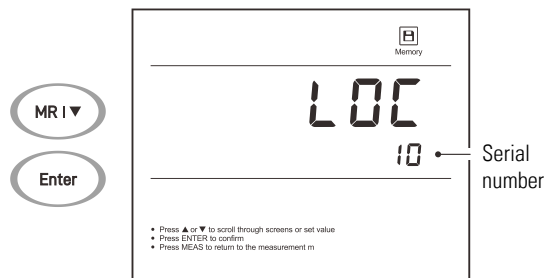

1.3 Press the ▼ key to view the stored data.

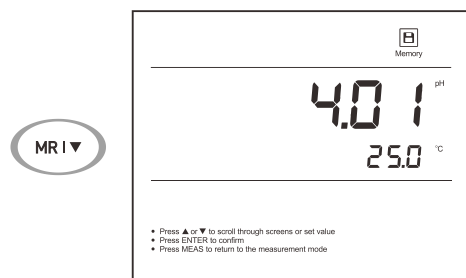

1.4 Press the ▼ key to view the next data set.

1.5 To exit the data log, press the **Meas** key.




If the meter does not store any reading, the display will show ---- only.

Clearing the Data Logs

If the memory is full, the meter will automatically show **FULL** when the **MI** key is pressed. To delete the data logs, please follow the steps below.

2.1 Press and hold the **MI** key to enter the setup menu.


2.2 Press the ▲ key until the meter shows **CLR/P-6**.

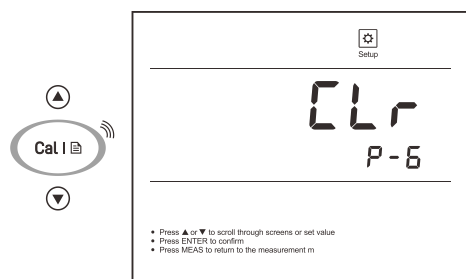

2.3 Press the **Enter** key, the meter shows **NO/CLR**.

2.4 Press the ▲ key to select the **YES/CLR**.

2.5 Press the **Enter** key to confirm.

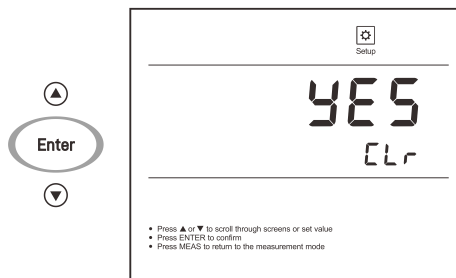


Communication

The 920E meter can transfer data to a computer or export data to Excel using DAS software. You can download the software from our official website. Before installation, make sure that Windows 10 is installed on your computer and that you have a USB-2303B data cable.

Receiving the Data

1. Connect the black port of the data cable to the meter and the blue port to the computer.

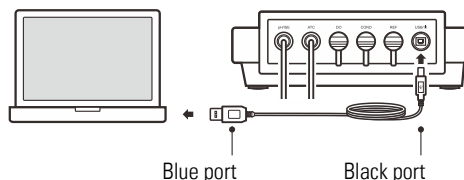

2. Double-click "**220_520_820_APP**". The system will automatically scan for an available port and display "Device port successfully recognized."
3. Click **OK** to start the application.
4. Click **Connect**. The screen displays "Connection established."
5. Click **OK** to continue.
6. Click **Receive**. The measurement values stored in the meter are automatically transferred to the computer.

Interval Recording

This function helps the user record measurement values at regular intervals over a specified time period. The recording interval can be set to 10/30/60 seconds, 10/30 minutes, or turned off.

1. Click the Interval Recording dropdown and select a parameter.
2. Click **Receive** to begin automatically transferring the measured values to the data sheet.

Creating an Excel File

1. Click **Save as Excel** to open the Save As dialog.
2. Enter a file name, then click **Save**.



Note: Once the software is closed, all received measurement values will be lost and cannot be recovered.

Electrode Maintenance

Cleaning the pH Electrode

Since pH electrode is susceptible to contamination, thoroughly clean as necessary after each use.

- **General Cleaning**
Rinse the pH electrode with distilled water and soak in 3M KCl solution.
- **Salt Deposits**
Dissolve the deposit by immersing the electrode in warm tap water. Rinse the electrode with distilled water and soak in 3M KCl solution.
- **Oil or Grease**
Place the electrode in the detergent or ethanol solution for 15 minutes. Rinse the electrode with distilled water and soak in 3M KCl solution.
- **Protein**
 - (1) Add 1% pepsin to 0.1M HCl solution.
 - (2) Place the electrode in above solution for 15 minutes.
 - (3) Rinse the electrode with distilled water and soak in 3M KCl solution.
- **Clogged Liquid Junction**
 - (1) Heat a diluted KCl solution to 60°C (140°F).
 - (2) Place the electrode into the heated solution for 10 minutes.
 - (3) Allow the electrode to cool in unheated KCl solution.

Reactivating the pH Electrode

If the pH-sensitive membrane has dried out, the electrode response will become sluggish. Immerse the electrode in a pH 4.01 buffer solution for about 30 minutes to rehydrate. If this fails, the electrode requires activation.

1. Soak the electrode in a 0.1M of HCl for 10 minutes.
2. Remove and rinse with distilled water, then place into a 0.1M of NaOH for 10 minutes.
3. Remove and rinse again, and soak in 3M KCl solution for at least 6 hours.

If these steps fail to restore the response, replace the electrode.

Storing the pH Electrode

- For best results, always soak the electrode in 3M KCl solution.
- DO NOT store the electrode in distilled or deionized water that will deplete the hydration layer of the pH-sensitive membrane and render the electrode useless.
- If you do not use the electrode for a period longer than 1 month, store the electrode in storage solution.

Cleaning the ORP Electrode

- Rinse the ORP electrode thoroughly with distilled water after use.

- In the corrosive chemicals, viscous solutions and solutions with heavy metals or proteins, take readings quickly and rinse electrode immediately.
- If the electrode response becomes sluggish, refer to the instructions below to clean the electrode.

- (1) **Inorganic Deposits**
Place the electrode in 0.1M HCl solution for 10 minutes. Rinse the electrode with distilled water and soak in 4M KCl solution.
- (2) **Oil or Grease**
Place the electrode in detergent such as dishwashing liquid for about 30 minutes. Rinse the electrode with distilled water and soak in 4M KCl solution.
- (3) If the platinum sensing element is severely contaminated, polish the platinum surface gently with an abrasive paper of 600 grid. Place the electrode in 0.1M HCl solution for 10 minutes. Remove and rinse with distilled water, then soak in 4M KCl solution for at least 6 hours.
If the electrode does not restore normal performance, replace the electrode.

Storing the ORP Electrode

If you do not use the electrode for long periods, store the electrode in 4M KCl solution.

Appendix

Optional Accessories

pH Electrodes

Order Code	Description
E201-BNC	For general applications
E202-BNC	For measuring semisolid and gel samples
P11	For measuring general and mildly corrosive liquids
P11-LiCl	For measuring non-aqueous liquids
P11-HA	For measuring high alkalinity liquids
P13	For measuring micro-volume samples
P15	For measuring low conductivity liquids
P16	For measuring samples with Tris buffer
P19	For piercing and measuring soft solids
P21	For measuring colloids

ORP Electrode

Order Code	Description
501	For the sample with strong redox potential
502	For the sample with weak redox potential
504	For the high temperature samples (< 80°C/176°F)

Temperature Probe

Order Code	Description
TP-10K	Range: 0 to 100°C (32 to 221°F), 1 m (3.3 ft.) cable

Solutions

Order Code	Description
PHCS-USA	pH 4.01, 7.00, 10.01 buffer solutions, 480 ml
PHCS-NIST	pH 4.01, 6.86, 9.18 buffer solutions, 480 ml
PHCS-ES	Electrode storage solution, 480 ml
PHCS-A	For removing acidic deposits, 480 ml
PHCS-B	For removing bacterial contaminants, 480 ml
PHCS-G	For removing oil and grease, 480 ml
PHCS-O	For removing organic contaminants, 480 ml
PHCS-P	For removing protein residues, 480 ml

Communication and Power Supply

Order Code	Description
USB-2303B	USB connector A to B, 1 m (3.3 ft.) cable
DCPA-5V	DC 5V power adapter, european standard plug

Preparation of Electrode Storage Solution

- Dissolve 24.6 grams of analytical grade potassium chloride (KCl) reagent in 100 ml distilled water.
- Add pH 4.01 standard buffer and adjust solution to pH 4.

Meter Specifications

Model	Bante 920E
pH	
Range	-2.00 to 20.00 pH
Resolution	0.01 pH
Accuracy	±0.01 pH
Calibration Points	1 to 3 points
pH Buffer Options	USA (pH 4.01, 7.00, 10.01) NIST (pH 4.01, 6.86, 9.18)

Auto-Buffer Recognition	Yes
Slope/Offset Display	Yes
Temperature Compensation	0 to 100°C (32 to 212°F), manual or automatic
mV	
Range	-1999 to 1999 mV
Resolution	0.1 (0.0 to 999.9), 1 (1000 to 1999) mV
Accuracy	±0.5 mV
Temperature	
Range	0 to 105°C (32 to 221°F)
Resolution	0.1°C (0.1°F)
Accuracy	±0.5°C (±0.9°F)
Calibration Point	1 point
Other Specifications	
Memory	100 data sets
Communication Interface	USB-B
Operating Temperature	0 to 50°C (32 to 122°F)
Storage Temperature	0 to 60°C (32 to 140°F)
Relative Humidity	< 80% (non-condensing)
Display	LCD, 125 × 100 mm (4.9 × 3.9 in.)
Power Requirements	DC 5V/400mA power adapter
Auto-Off	30 minutes after last key pressed
Dimensions	210 (L) × 188 (W) × 60 (H) mm, (8.2 × 7.4 × 2.3 in.)
Weight	1.5 kg (3.3 lb)

Troubleshooting

Fault	Cause and Corrective Action
Screen shows - - - -	Electrode dried out. Soak the pH electrode in 3M KCl solution for about 30 minutes.
	Measurement exceeded the maximum range. Check the electrode and sample.
Drifting erratic readings	Check whether electrode is contaminated, clogged or broken.
Screen shows E r r	pH buffer problem. Use freshly prepared buffer solutions to calibrate the meter.

Disposal

This product is required to comply with the European Union's Waste Electrical and Electronic Equipment (WEEE) Directive 2002/96/EC and may not be disposed of in domestic waste. Please dispose of product in accordance with local regulations at the collecting point specified for electrical and electronic equipment.



Warranty

The warranty period for meter is one year from the date of shipment. Above warranty does not cover the electrode and pH buffer solutions.

Out of warranty products will be repaired on a charged basis.

The warranty on your meter shall not apply to defects resulting from:

- Improper or inadequate maintenance by customer
- Unauthorized modification or misuse
- Operation outside of the environment specifications of the products

For more information, please contact the supplier.



Office: 4715 Castlewood St., Sugar land, TX 77479, USA

Tel: (+1) 346-762-7358

E-mail: banteinstruments@yahoo.com

Factory: F3, Building 2, No.2185, Laifang Rd., Shanghai 201615, China

Tel: (+86) 21-6404-1598

E-mail: banteinstrument@hotmail.com

